

12/14 Morse Taper

Collar Design
prevents stem sinking

Shark-Fin Design
pressurizes cement during
stem implantation

CoCrMo Forging Stem
high fatigue strength

Tri-Wedge Design
facilitates proximal
load transfer

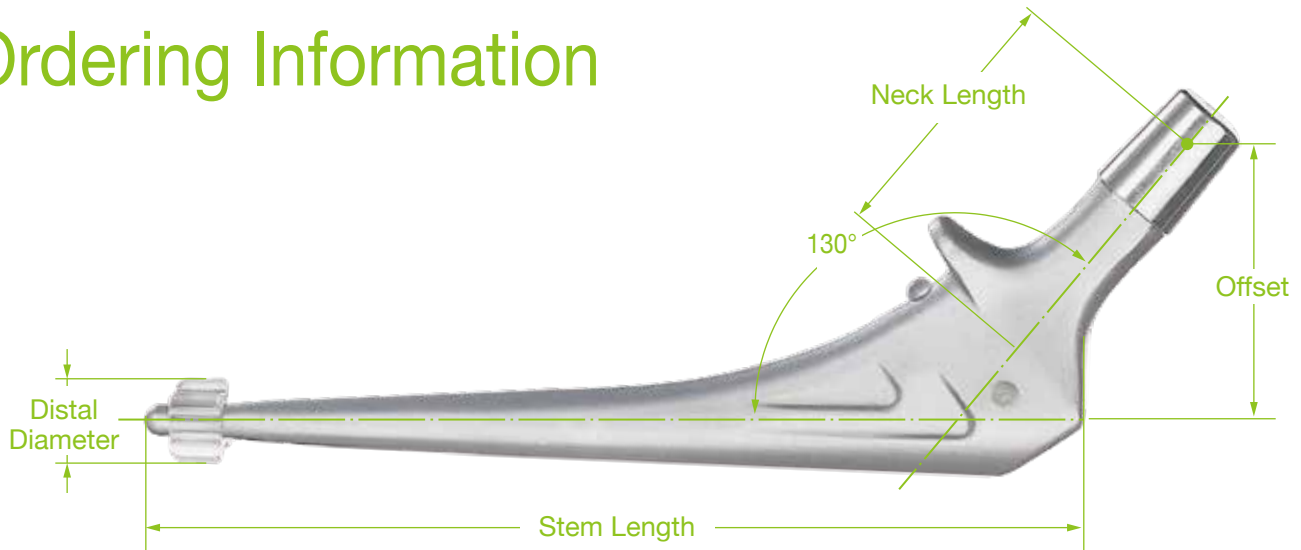
2 mm PMMA Spacer
for even cement thickness

U2[™] Hip System Cemented Stem

Distal Centralizer
equalizes cement mantle and
ensures correct stem alignment


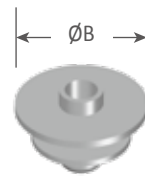


Ordering Information



Catalog Number	Stem size	Distal Diameter (mm)	Stem Length (mm)	Cement Mantle (mm)	Offset (mm)	Neck Length (mm)
1104-5223	#2	Ø10	126	1.0-2.0	36	30
1104-5234	#3	Ø11	132	1.5-2.0	39	34
1104-5245	#4	Ø12	138	1.5-2.0	39	34
1104-5256	#5	Ø13	144	1.5-2.0	43	38
1104-5267	#6	Ø14	151	1.5-2.0	45	41

Accessories

Catalog Number	Description	(mm)		NEW	Catalog Number	Description	(mm)	
		ØA	ØB				SIZE	ØB
 Cement Restrictor	1905-7007	7	13	 Cement Restrictor, I-Type	1907-1008	8~10	12	
	1905-7008	8	14		1907-1010	10~12	14	
	1905-7009	9	15		1907-1012	12~14	16	
	1905-7010	10	16		1907-1014	14~16	18	
	1905-7011	11	17		1907-1016	16~18	20	
	1905-7012	12	18		1907-1018	18~20	22	
	1905-7013	13	19					
	1905-7014	14	20					
1905-7015	15	21						
1905-7016	16	22						



Each Step We Care



Contact Us

■ Taiwan

United Orthopedic Corporation

12F, No. 80, Sec. 1, Chenggong Rd, Yonghe Dist.,
New Taipei City 23452, Taiwan

Tel: +886 2 2929 4567 Fax: +886 2 2922 4567

■ Europe

UOC (Suisse) SA

Avenue Général Guisan 60A, 1009 Pully, Switzerland.

Tel: +41 (0) 21 634 70 70 Fax: +41 (0) 21 634 70 76

UOC (France)

21, Rue de la Ravinelle, 54000 Nancy, France.

Tel: +33 (0) 3 83 23 39 72 Fax: +33 (0) 3 83 23 39 10

■ China

United Medical Instrument Co. Ltd.

Room 402, Building 7, No.697, Ling ShiRoad,
Shanghai, China, 200072

Tel: +86 21 3661 8055 Fax: +86 21 6628 3031

■ EU Representative

mdi Europa GmbH

Langenhagener Strasse 71,
30855 Langenhagen, Germany

Tel: +49 511 3908 9530 Fax: +49 511 3908 9539

■ Japan

United Biomech Japan Inc.

Ginyo Bldg 4th floor, 2-9-40 Kitasaiwai, Nishi-Ku
Yokohama City, Kanagawa, 220-0004, Japan

Tel: +81 45 620 0741 Fax: +81 45 620 0742

■ USA

UOC USA Inc.

20 Fairbanks, Suite 173, Irvine, CA 92618

Tel: +1 949-328-3366 Fax: +1 949-328-3367